



Frontline Staff Quick Reference Guide

*Disinfection Tool Kit
for Pre-Operative and PACU Areas
Using Hydrogen Peroxide Wipes*



The AORN Seal of Recognition was awarded to the Clorox Healthcare Disinfection Tool Kit for the Operating Room, Pre-Operative and PACU Areas on November 26, 2019 and does not imply that AORN approves or endorses any product or service mentioned in these materials. The AORN Seal of Recognition program does not include any CE credit for programs.





1. Stick to one cleaning path: clockwise or counterclockwise.
2. Use the Zone Approach to clean.
3. Clean from high to low.
4. Wipe from clean to dirty.

7 Steps for Success



PREPARE

Gather supplies, perform hand hygiene and put on PPE.



INSPECT

Inspect the surface to see if there is any visible soil or foreign material.



REMOVE SOIL (CLEAN)

If there is soil or foreign material on the surface, you will need to clean first with one wipe.



WIPE SURFACE (DISINFECT)

Using a fresh wipe, apply to the surface in an S-pattern, moving from clean to dirty and overlapping slightly on each sweep of the surface. Do not go over the same area twice with the same wipe. Use more than one wipe on large surfaces if necessary.



WAIT (CONTACT TIME)

Ensure the surface remains visibly wet for the full contact time listed on the product label.



DISCARD

Discard the used wipe when you have finished wiping the surfaces, or if the wipe becomes soiled or dry.



POLISH

Wipe surface with a clean, dry or damp cloth to prevent residue build-up and protect surfaces.

Hydrogen Peroxide Cleaner Disinfectant Wipes



Products & Safety

Product Overview

Clorox Healthcare® Hydrogen Peroxide Cleaner Disinfectant Wipes

CONTACT TIME:*



DIRECTIONS FOR USE



1. Don PPE.
REMOVE gross,
visible soil.



2. APPLY directly
onto surface
evenly until
thoroughly wet.



3. To DISINFECT,
keep surface
wet for 1 minute.



4. POLISH
surfaces to
remove residue
if visible.

At a Glance Hydrogen Peroxide Cleaner Disinfectant Wipes

Active Ingredient	Hydrogen Peroxide
Format	Pre-moistened wipes
Dilution	Ready to use (no mixing required)
Contact time*	1 minute
Kills <i>C. difficile</i> spores?	No
Shelf Life	Expiration date listed on each package
PPE Required	Follow your facility's protocol

Refer to Safety Data Sheets in your EVS Manager's Toolkit for complete information.

* Use as directed, refer to the EPA master label for complete list of organisms & contact times.

Use products as directed. Clorox Healthcare recommendations.
Your facility's protocol may vary.
© 2019 Clorox Professional Products Company.

NI-48554 Updated November 2019





About the product

- This is a one-step cleaner disinfectant product with hydrogen peroxide as the active ingredient.
- The formula also contains detergents and other ingredients to speed up the disinfecting power of hydrogen peroxide and help improve the product's cleaning power.
- These detergents are found in common cleaning products.
- While hydrogen peroxide itself has no smell, the other ingredients make the product smell like a common cleaner.



Residue and Surface Compatibility

- Hydrogen peroxide wipes are fast-acting and effective cleaners and disinfectants.
- Hydrogen peroxide disinfectants can be used facility-wide on a wide variety of surfaces including stainless steel and plastics.
- Remember that all disinfectants can leave some residue on surfaces.
- If you see some residue left behind, you can wipe it away with a clean, dry or damp cloth.


Pre-Operative and PACU:

High-Touch Surfaces



Protocols

Hydrogen Peroxide Wipes

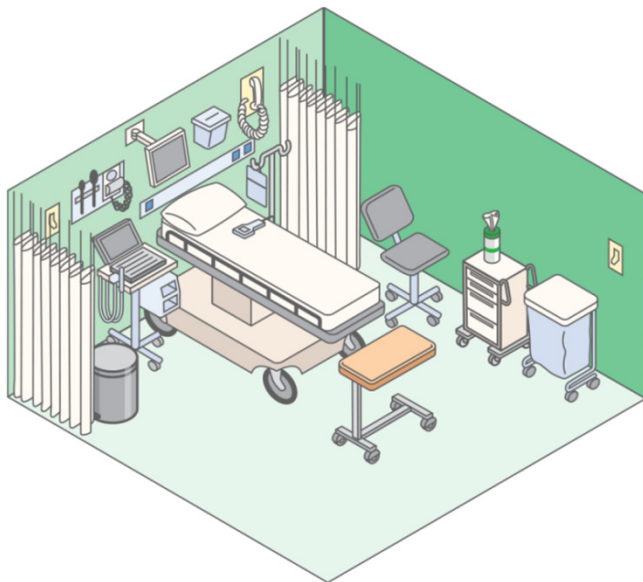

CLEAN AND DISINFECT HIGH-TOUCH SURFACES WHICH MAY INCLUDE:



Bedside table



Glucometer



Patient warming device



Call button



Infusion pump + IV pole



Phone



Cardiac monitor



Light switches



Trash can



Computer on wheels



Overbed table



Wall panel behind bed



Door handles/
push plates



Patient bed



Vital signs monitor



Use products as directed. Clorox Healthcare recommendations.
Your facility's protocol may vary.
© 2019 Clorox Professional Products Company.


NI-48554 Updated November 2019

Pre-Operative and PACU Areas:

Between Patients & Nightly



Hydrogen Peroxide Wipes

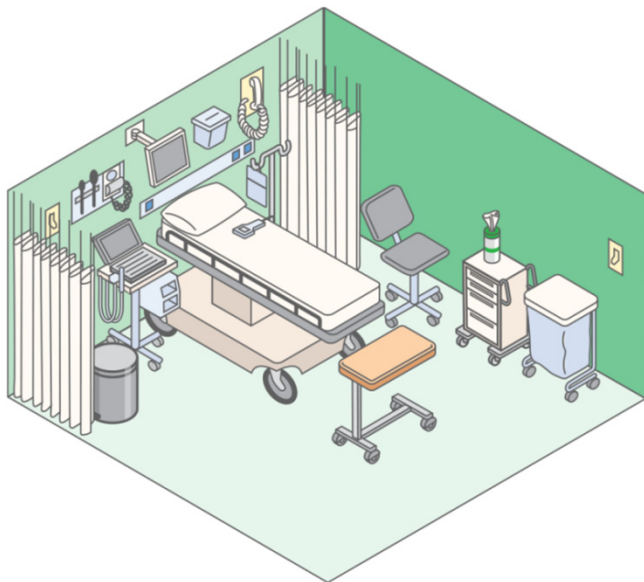

PRE-OPERATIVE/PACU: BETWEEN PATIENTS & NIGHTLY

- Wear the appropriate PPE including gloves per facility protocol.
- Clean to remove gross soil before disinfecting.
- Clean from high to low and from clean to dirty.
- Clean and disinfect monitors, bed, overbed table, TV remote, call light and equipment, including all high-touch surfaces.
- Use more than one wipe as needed for larger surfaces.
- Keep the surface visibly wet for the full contact time .
- Mop floors using disinfectant between patients only if soiled.
- Mop floors with disinfectant nightly.